
DERMASUCTION™

Pore Cleaning Device

IMPORTANT SAFETY INSTRUCTIONS, READ AND FOLLOW ALL WARNINGS AND INSTRUCTIONS BEFORE USING THIS PRODUCT. INJURY CAN RESULT FROM IMPROPER USE. SAVE THESE INSTRUCTIONS.

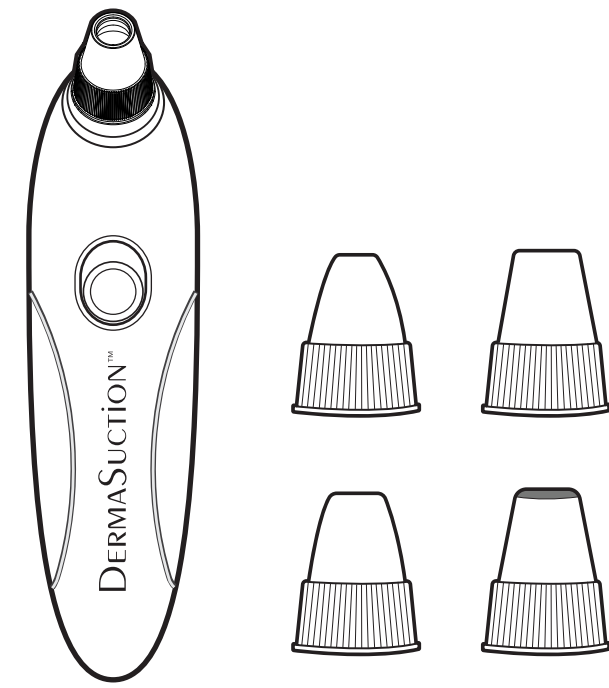PRODUCT MAY VARY FROM IMAGES SHOWN

TABLE OF CONTENTS

Warnings	1
Product Components	2
Battery Installation	3
How To Change The Suction Heads	3
How To Use	4
How To Clean	5

This AC wire carrier current device of unintentional radiator is a Category II equipment compliant with Industry Canada's ICES-006.



WARNING

- Before use, talk to your doctor if you have any medical issues or concerns.
- This product is not for use by children.
- Necklaces, earrings and other types of jewelry should be removed before using this product.
- This product is not for use on areas of the body that have undergone any sort of surgery or have any scars.
- If you experience pain while using this product, discontinue use immediately.
- This product is not water proof and should not be used in the shower or any situation that involves water.
- Do not store this product in the bathroom or a room that has high moisture content.
- Do not operate in the presence of explosives and/or other flammable fumes or liquids.
- Do not store or use this product near sources of heat or fire.


WARNING

KEEP BATTERIES AWAY FROM CHILDREN. HARMFUL IF SWALLOWED.


PRODUCT COMPONENTS

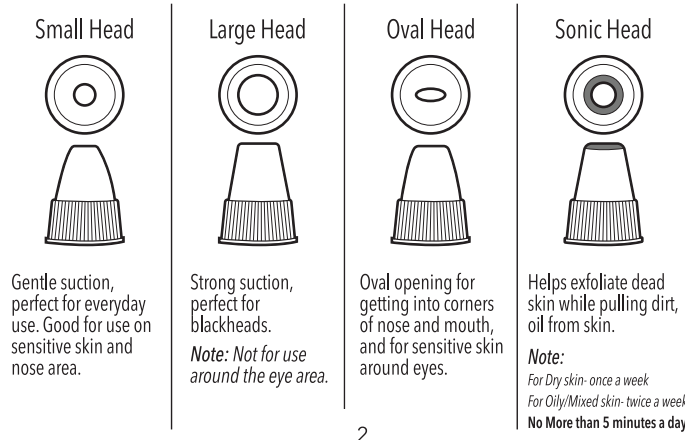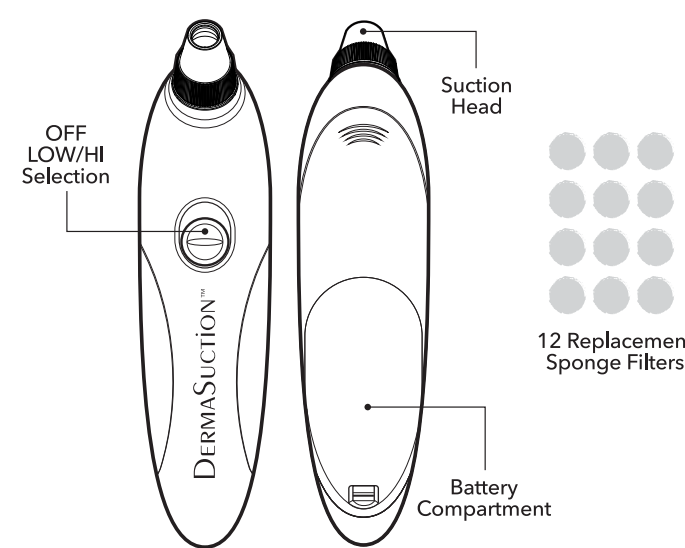
BATTERY INSTALLATION

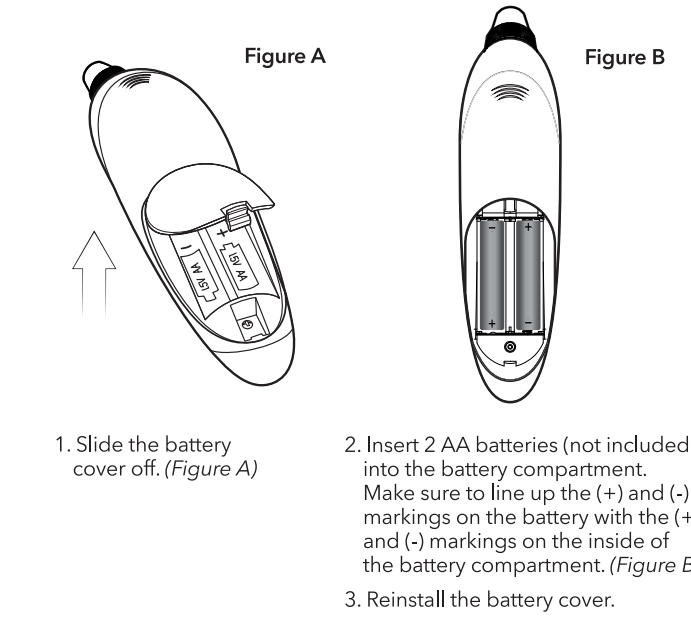HOW TO CHANGE THE SUCTION HEADS

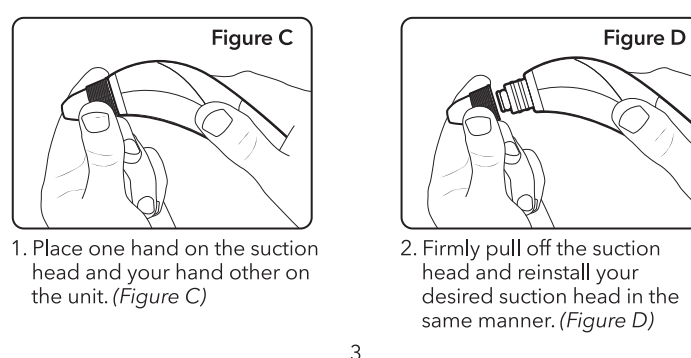
HOW TO USE

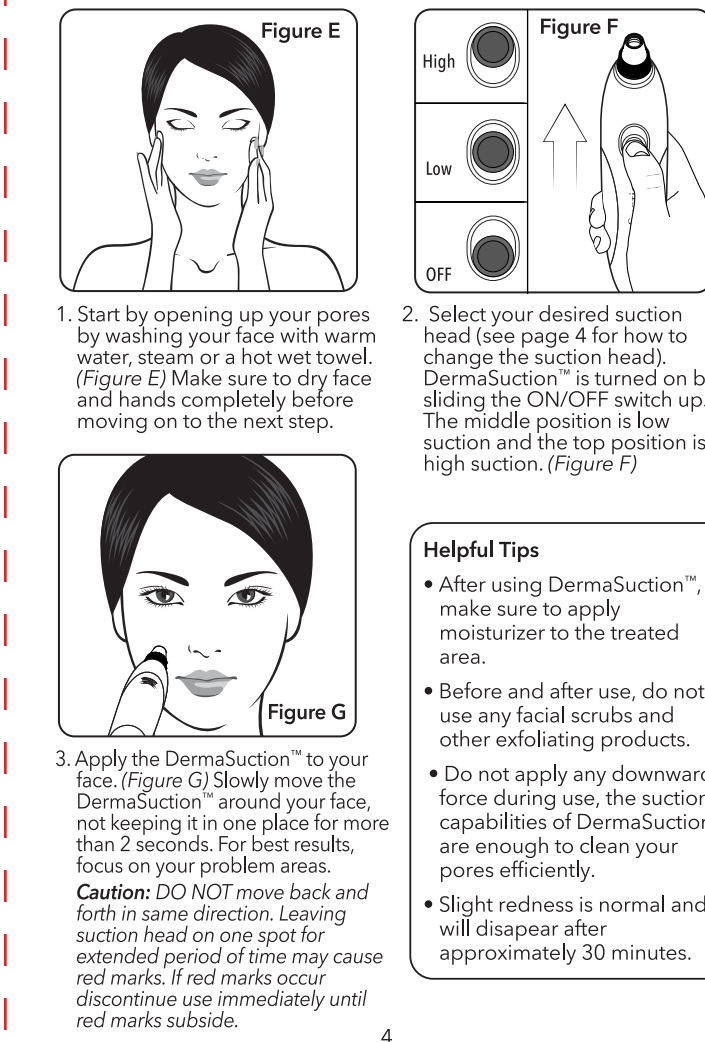- Helpful Tips**
- After using DermaSuction™, make sure to apply moisturizer to the treated area.
 - Before and after use, do not use any facial scrubs and other exfoliating products.
 - Do not apply any downward force during use, the suction capabilities of DermaSuction™ are enough to clean your pores efficiently.
 - Slight redness is normal and will disappear after approximately 30 minutes.

HOW TO CLEAN

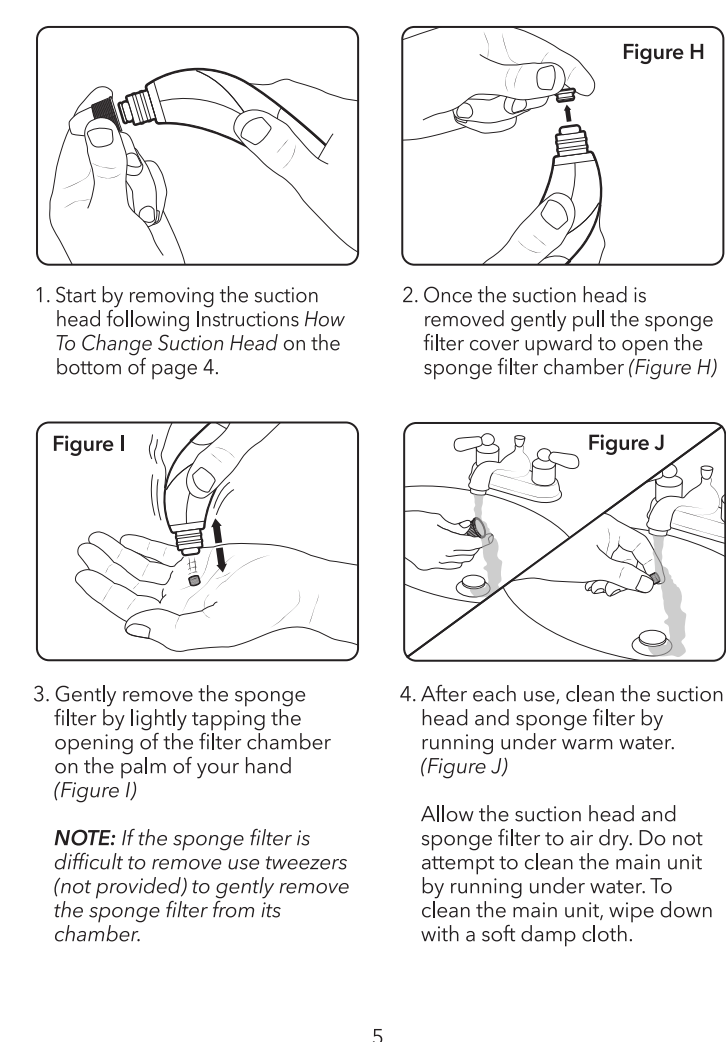Purchase additional Sponge Filters here:
www.clickit.ca

LIMITATION OF LIABILITY
Liability limited to the purchase price of this product. Supertek Canada Inc. shall not be liable for any incidental or consequential damages or any implied warranty on this product.

Distributed by
Supertek Canada Inc.
Mont-Royal, Qc H4T 1X2
Customer Service:
1-800-304-8354
www.clickit.ca

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