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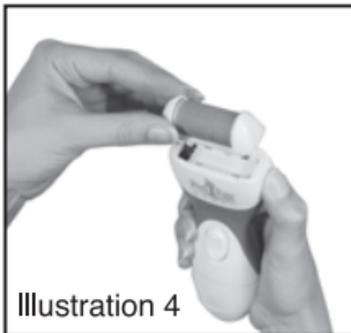
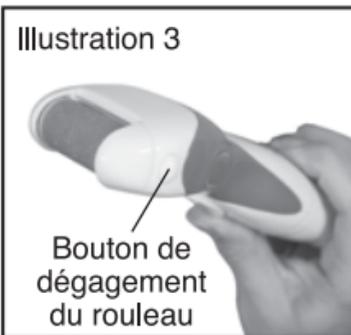
CAUTION:

- Intended for use on callused feet only. Do not use on uncallused, smooth skin. Injury may result.
- Avoid use on irritated, inflamed, or infected skin, or on open wounds. Irritation or laceration may result.
- Consult a physician before use if any skin conditions exist.
- Avoid use if you are diabetic or have poor blood circulation. Consult a physician before use.
- PedEgg Power is not recommended to be shared or used by more than one person. Infection may result.
- Not recommended for use by children.

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LIMITATION OF LIABILITY


Liability limited to the purchase price of this product. Supertek Canada Inc. shall not be liable for any incidental or consequential damages or any implied warranty on this product.

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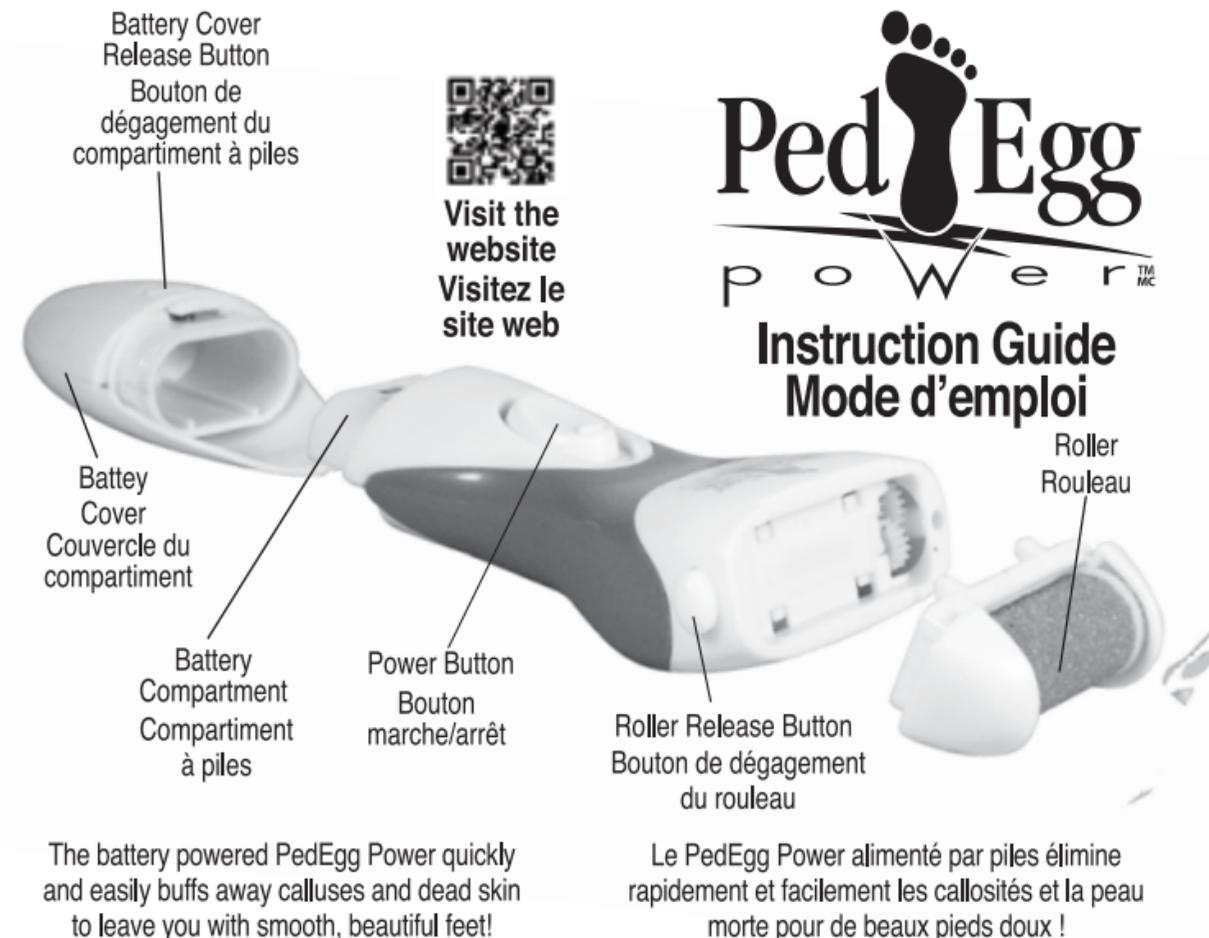
⚠ WARNING:

Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (nickel-cadmium) batteries.

Before First Use:

1. Remove battery cover by pressing release button on the underside of the PedEgg Power then pulling the battery cover downwards. (Figure 1)
2. Insert 2 AA batteries.
3. Close battery cover by pressing battery cover up until release button locks into place.

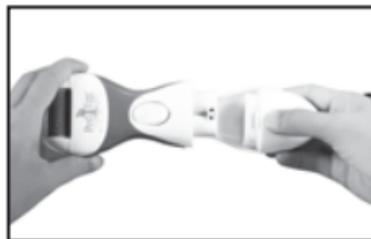**Instructions For Use:**

1. Turn on the PedEgg Power by pressing the power button up.
2. Holding the PedEgg Power securely, gently press the roller to the thickest area of dead skin you want to remove and slowly move outwards towards healthy skin. (Figure 2)
3. Check intermittently for desired results. Once new, healthy skin is revealed or callus is removed, stop use and switch the PedEgg Power off.
4. For additional smoothing, switch from the blue

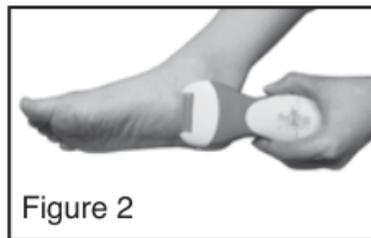
Figure 2

NOTE: The PedEgg Power includes an auto off feature which stops the roller if you are pressing too hard.

Nano Abrasion Roller to the purple Smoothing Roller (see *To Remove Roller below*) and use gentle pressure to smooth the areas you buffed.

NOTE: Stop immediately if skin becomes irritated or if any pain is felt. Overuse may cause skin irritation and soreness.



To Remove Roller:

1. Press the button on the side of the PedEgg Power near the roller. (Figure 3)
2. Pull upwards on the roller to remove it. (Figure 4)
3. Snap new roller into place by pressing it down onto the PedEgg Power.

Cleaning Instructions:

1. With PedEgg Power switched off, remove roller following instructions *To Remove Roller*.
2. Rinse roller under running water.
3. The body of the unit can be wiped clean with a dry cloth.
4. Allow roller to dry and insert it back into PedEgg Power body as instructed in *To Remove Roller*.

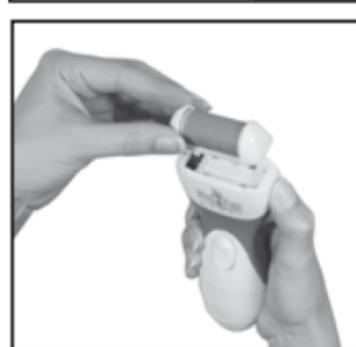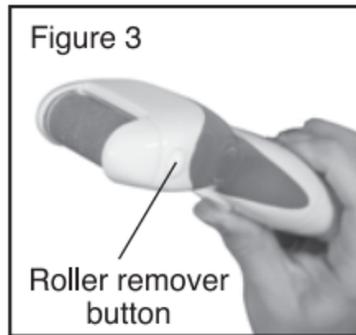NOTE: For best performance, clean PedEgg Power after each use.

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